



TR 4000-W Atlas

Battery operated shower trolley

Hygiene Equipment
Designed for People

TR 4000-W Bariatric shower trolley

TR 4000-W Multi tasking bariatric battery operated shower trolley is unique with it's double columns, allowing the trolley stretcher to overlap the bed for a safe and easy patient transfer.

The trolley drains in both ends, allowing the person to be positioned in either direction, an invaluable feature for the caregiver.

The shower trolley is available with an onboard battery solution, or with removable rechargeable battery.

Central braking and straight steering castors is standard, it will make the trolley easier to maneuver.

The height adjustable shower trolley is a great solution to prevent back injury and minimize static load. It offers a safe and ergonomic working environment for the caregiver.



Registered at Medical Products Agency, Sweden as a MDR Class 1 product according to MDR 2017/745.

TR reserves the right to make design changes and/or improvements without notice.





TR 4000-W Atlas

Battery operated shower trolley

Hygiene Equipment
Designed for People

- Quiet and smooth height/tilt adjustable motion
- Designed for easy cleaning and ease in mobility
- Mattress easily removable
- Foldable side rails for efficient and safe patient transfer and transport.
- Dual drainage offers easy changing of head end.
- The side rails can be pressed together for traveling through narrow doorways.
- The stretcher can easily be tilted in both directions to a trendelenburg position. An audible tone will sound when the stretcher is level.
- The chassis and side frames are made of powder coated stainless steel.
- The stretcher has rubber bumper wheels on all four corners to protect doors and doorways.
- Proven choice of durable high quality components guarantees a low cost of ownership.
- Both versions come with a magcode charger, mattress, head pillow, 4 castors and 1 m flexible drainage hose with hose holder, as standard.

Specifications

Length	2205 mm
Width	1067 mm
Height	535 - 835 mm
Floor clearance	160 mm
Working load limit	450 kg
Castors	4 central braked castors, 150 mm as standard.
Control	Hand control for height and tilt adjustability.
Battery	Onboard 24 V, 7 Ah / Removable Li-Ion 24 V, 4 Ah.
Protection class	Complete unit IPX5 (Hand control and removable battery IPX6)

Models:

TR 4000-W 2205 x 1067 mm


Options/Accessories:

Backrest wedge pillow 15° or 30°

Fasteners/straps for attachment on mattresses to siderails

Longer drainage hose 2.5 m

Removable rechargeable battery solution.

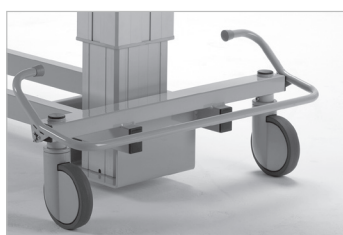

Central breaking/steering



Head end rail folds for extra length

